

Beauty Equipment



*We help You
to look perfect...*

Aqua-jet peeling - is an innovative technology of highly effective and pain-free skin rejuvenation. The treatment effect is based on non-invasive infusion of compressed oxygen and medical solution such as saline or special serums with necessary vitamins, microelements, collagen or hyaluronic acid)

Treatment Description

Aqua-jet peeling treatment combines 100% oxygen and sterile saline to gently cleanse and hydrate the skin. This solution under the high pressure comes out of the nozzle and treat the necessary skin area. Oxygen penetrates into the deep dermis and provides cell oxygenation treatment. Special serum (if it is used during the treatment) penetrates into dermis and nourishes it. Compressed oxygen and medical solution provides massage treatment together with lymph drainage and lifting effect. It exfoliates the top layers of dry dead skin cells and clean out the skin surface.


While choosing the distance between the nozzle tip and the skin, it is possible to regulate the intensity level and the treatment result: from gentle massage effect to the deep skin peeling (much like a microdermabrasion treatment). One treatment movement provides stimulation effect, improves blood circulation and micro circulation, causing a characteristic blush.

Several movements, on the contrary, cause the outflow of blood from the capillary vessels located close to the surface of the skin (which is excellent remedy to fight the first signs of Rosacea the repeated close to skin movements

Aqua-jet peeling AV – 1100



Aqua-jet peeling AV – 3000

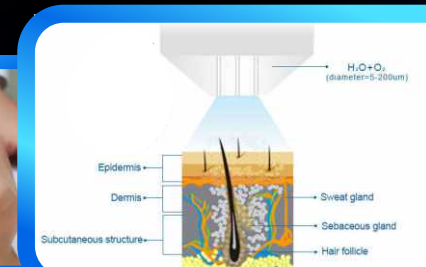

Advantages:

- mesotherapy without needles and injections
- clean technology
- 100% hygienic contactless treatment
- absence of chemicals and harmful substances
- massage effect, providing blood circulation improvement
- pain-free and 100% comfort
- no downtime
- immediate visible results



Multifunctional beauty machine N-02



Function:

- needle-free mesotherapy
- diamond dermabrasion

Multifunctional beauty



Function:

- ultrasonic scrubber
- diamond microdermabrasion
- needle-free mesotherapy (electroporation)

Multifunctional beauty



Functions:

- ultrasonic scrubber
- diamond microdermabrasion
- hot-cold therapy
- photon therapy
- needle-free mesotherapy (electroporation)

Multifunctional beauty



Functions:

- ultrasonic scrubber
- diamond microdermabrasion
- hot-cold therapy
- photon therapy
- needle-free mesotherapy (electroporation)
- microcurrent therapy

Multifunctional beauty machine F-4



Functions:

- ultrasonic scrubber
- needle-free mesotherapy
- microcurrent therapy


Parameters:

- oscillation frequency of ultrasonic blade: 27-28 kHz
- electroporation: 2,5-5 kHz, 36W, 10-30mA
- hot-cold hammer:
 - hot electrode: from 40°C up to 55°C
 - cold electrode: from 5°C up to 20°C
- spot size: 5,5 cm
- output microcurrent: 10-90 mA
- microcurrent modes: 4 modes
 - output microcurrent:
 - bi polar pulses with harmonic modulation
 - diamond dermabrasion: 2 nozzles,
 - 9 nozzle tips with diamond coating
 - vacuum suction in diamond dermabrasion: 60-70 cm hg.
 - photon- 7 lights

Assemble your own multifunctiona beauty machine



Needle-free mesotherapy machine T-12

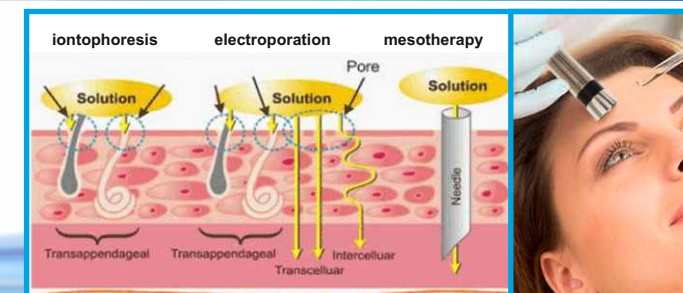

The main advantage of needle-free mesotherapy for patients is the absence of painful sensations without loss of therapeutic effect. In addition, this type of treatment is widely used by cosmetologists in beauty salons because of its multifunctionality. Only the one beauty machine is needed to eliminate a variety of different skin blemishes of the face and body. This makes it must-have of any cosmetologist.

Indications:

- aging skin
- excessive sensitivity and dryness of the skin
- the tendency to inflammation
- cellulite
- pigment spots
- hair loss
- acne rash
- postnatal skin defects



The operation principle is based on microcurrents that penetrate into the deep dermis and favorably effect on it. It facilitates facial and décolleté area skin health, production of proteins, lipids and enzymes, the restoration of the skin's ability to regenerate. T-15 beauty machine also influences mimic muscles and has the effect of its stimulation. During the treatment, the muscles are continuously compressed and decompressed, thereby affecting the vascular system.

Parameters:

- Output current: 20-700 mA
- Pulse frequency: 1-400 Hz
- Intensity regulation: from 1 to 15
- 4 operation modes:
 - P1 – reprogramming of mimic muscle
 - P2 – micro iontophoresis (ionic mesotherapy)
 - P3 – lifting of face and neck (it can be used for body areas as well)

Vacuum therapy beauty machine T-01



When using T-01 vacuum therapy beauty machine, the skin health improves significantly, its elasticity is tightened and restored, lymph drainage is increased and blood circulation is improved. This method is also used to improve cellulite.

Parameters:

- Frequency: 50 Hz - 60 Hz
- Vacuum suction: 60-70 mm Hg.
- Set of delivery: 3 roller nozzles
- Operation modes: pulse and continuous

Ultrasonic scrubber beauty machine T-02



The operation principle is based on the influence of ultrasonic that helps to gently peel off the keratinized cells and dirt, to improve pH balance metabolic process in the skin tissue.

Functions:

- ultrasonic peeling
- phonophoresis



Electrocoagulation is one of the most popular and safe method to remove neoplasms that appear on the skin surface. Such way of treatment allows removing neoplasm and coagulation blood vessels at the same time. Electrocoagulation is one of the delicate methods, as, unlike surgical intervention, after it the risk of formation of a coarse rumen on the skin is minimal. There are no scars on the skin surface after such kind of removal. Such aesthetic treatment is used to remove neoplasm or papilloma on sensitive areas of ciliary margin, canthus etc.

Galvanic therapy beauty machine T-04



- 2 pcs set of nozzles (roll, round)
- output voltage: 1-30 B
- output current electroplating: 0,1 - 6,5 mA

Hot-and-cold therapy beauty machine T-05




- temperature regulation
- temperature of hot electrode: 40 C up to 55 C
- температура холодного электрода: от 5 C до 20 C

Diamond dermabrasion beauty machine T-06



Diamond dermabrasion is one of the most effective and safe method of mechanical peeling. T-06 beauty machine combines painlessness and effectiveness. It helps to treat different skin blemishes. The treatment is based on exfoliation of keratinized skin cells with the help of special abrasives - diamond crystals.

The unique system of T-06 facilitates such aesthetic skin blemishes as fine lines and wrinkles, acne, skin spots. It prepares the skin to other types of skin treatments.

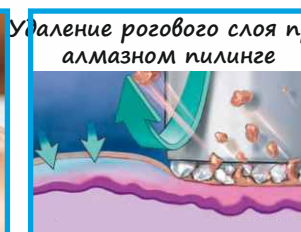

Set of delivery:

- 2 nozzles,
- 9 nozzle tips with diamond coating



Brushing therapy beauty machine T-09



The purpose of the rotary brush is to peel off the keratinized cells. It is used for activating massage and skin cleansing, removing of keratinized cells from the skin surface. Brushing massage helps to improve blood circulation.

- big brush nozzle (1 pc.)
- small brush nozzle (1 pc.)
- big pumice nozzle (1 pc.)
- small pumice nozzle (1 pc.)



High Frequency Beauty Machine T-10



- Pulse repetition frequency: 100 ± 10 Hz
- (Min): 8 kV peak value
- output voltage (Max): 25 kV



Ultrasonic therapy beauty machine T-11



- ultrasonic intensity: 1; 3 MHz
- equipment:
 - flat nozzle (to face treatments, neck, neckline)
 - conical (for the eye area, lips, nose)

Electrocoagulation and electroepilation beauty machine T-18



- operating frequency: 2700 Hz and 1.65 MHz
- decrease / increase of current mode



RF lifting beauty machine T-14



- 2 RF multipolar nozzles (10 mm, 32 mm)
- power: 15-35 J/cm²
- pulse type: DC/AC
- frequency: 3 MHz



RF lifting beauty machine PX-3000



- intelligent temperature control
- LCD touch screen 7"
- 3 multipolar nozzles RF (10mm, 32mm, 50mm)
- frequency RF: 3-5 MHz

Radiofrequency lifting

Radiofrequency lifting (radiolifting, RF-lifting) is the innovative technology for effective rejuvenation and non-surgical skin tightening of face and body, which is an excellent alternative to surgical lifting.

This treatment stimulates the activity of fibroblasts, which are responsible for collagen production, molecules of glycosaminoglycan (in particular, hyaluronic acid) and elastin. As the result, collagen fibers become tight and elastic, and numerous metabolic processes and fibroblasts are activated. Under the influence of heat, young elastin is produced. This process leads to the formation of a new healthy frame, which effectively supports the skin and returns beautiful and youthful appearance.

The radio frequency face lift is used for:

- ptosis (sagged skin)
- aging skin
- crow's feet
- dull skin
- face oval correction double chin
- dark under eyes circles
- skin puffiness
- photo aging
- deep fine lines and wrinkles
- acne scars

RF-lifting of body areas (belly, hips, buttocks) has the following indications:

- cellulite (fibrous and edematous)
- loss of skin elasticity, sagged skin
- stretch marks
- excess fat deposits
- photoaging

PX-3000 beauty machine has an innovative system of temperature control. It helps to achieve maximum effect within safety limits. Temperature control indicator built in front panel of the machine helps to avoid skin burning and overheating. It makes the treatment safe and comfortable

Multifunctional beauty Multifunctional beauty Multifunctional beauty



- Functions:**
- ultrasonic scrubber
 - diamond dermabrasion
 - ultrasonic therapy (1,3 Hz)



- Functions:**
- ultrasonic scrubber
 - diamond dermabrasion
 - ultrasonic therapy (1,3 Hz)



- Functions:**
- ultrasonic scrubber
 - diamond dermabrasion
 - ultrasonic therapy (1,3 Hz)
 - hot-cold therapy

Needle-free mesotherapy beauty machine S-02 (2 & 1)



- Functions:**
- ultrasonic scrubber
 - mesotherapy

Ultrasonic therapy beauty machine S-03 (2 & 1)



- Functions:**
- ultrasonic scrubber
 - ultrasonic therapy (1,3 Hz)

Multifunctional beauty



- Functions:**
- ultrasonic scrubber
 - needle-free mesotherapy
 - microcurrent therapy (face and body)

Multifunctional beauty



- Functions:**
- ultrasonic scrubber
 - needle-free mesotherapy
 - microcurrent therapy
 - diamond microdermabrasion

Neodymium laser TL-500



Functions:


- carbonic peeling
- permanent make-up removal
- tattoo removal




TL-500 provides carbonic rejuvenation (carbonic peeling) of facial skin. Laser carbonic peeling combines the effect of laser and rejuvenation treatment. The basis of treatment is a special gel mask. It consists of nanoparticles of carbon dioxide, and therefore the Nanogel is called carbon. The laser selectively acts on the pigment by using a very high power pulse, without affecting nearby tissues. Color pigment tattoo is broken and successively removed by the immune system. Therefore, there are no negative effects after the treatment (cicatrices, scars, white spots).

TL-500 emits high energy that instantly destroys the pigment of tissues at the place of treatment. High energy acts instantaneously, radiating laser waves of certain length and penetrating the epidermis and falling into the pigment of the tissue. Then these pigments expand at maximum speed and are destroyed. Part of the pigments are divided into granules, which are naturally excreted through the lymphatic system. The pigment gradually grinds and disappears.

Laser

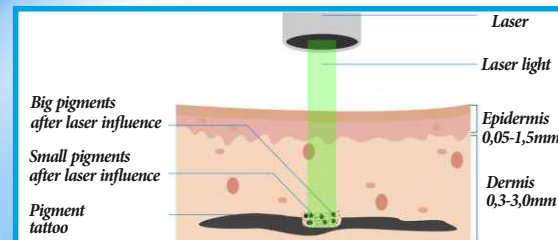

Tattoo removal process

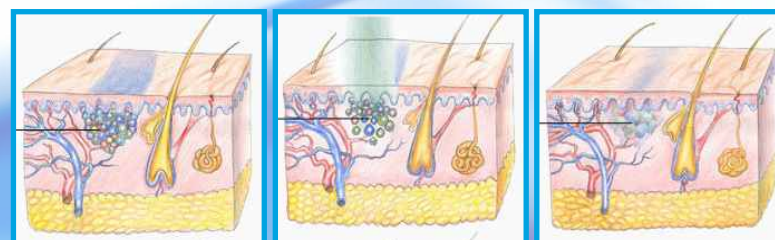


Pigment tattoos

laser light destroys molecules of pigment – destroying of skin cells

fragments of pigment – tattoo is decolorized

Tattoo removal stages

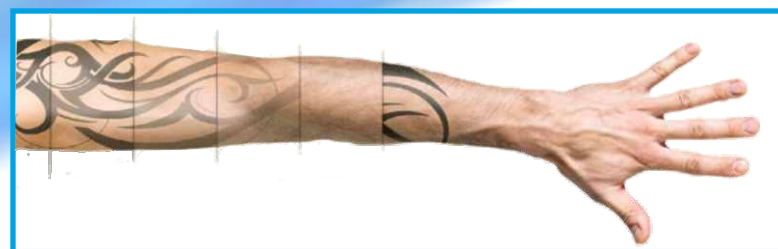

Parameters:

- Laser type: -Nd:YAG Q-Switched

Laser

- Wavelength: 1064nm, 532nm, 1320 nm (carbonic peeling)
- Energy: (max) 1000mJ
- Frequency: 1-10Hz, adjustable
- Power: 300-500W
- Pulse time: <10ns
- Voltage: 110-220V, 50-60Hz
- Spot diameter: 2-5mm
- Cooling system: Wind cooling and water cooling

ALVI
Frague

Diode laser ALVI DL - 6000

Esthetically Hygienically Permanently



Diode laser Alvi DL- 6000 is a leader in the field of laser systems for removal of unwanted hair.

Currently, diode laser technique remains the most effective, as it allows working on tanned skin, remove any hair thickness, color and density except vellus, very light and gray hair.



Functions:

- Hair removal

Parameters:

- Wavelength: 808 nm:
- Spot size: 15 * 15 mm
- Pulse Width: 10-400 ms
- Energy Density 100 J/cm²
- Pulse Frequency: 1-10 Hz
- Cooling Mode: cooling sapphire tip, air and water cooling system, semiconductor
- Cooling Temperature: 2C ~ -10 C
- Power output: 600 W
- Operating mode: HR, SHR
- Lamp life cycle: 10 000 000 laser shots




ALVI
Frague

Photoepilator IP-200

It combines three of the most popular photo technologies. The machine is equipped with two types of handpieces HR and SR, manufactured in Germany with the resource of 300,000 pulses for each one. The beauty machine is used for the removal of unwanted hair of all types, skin rejuvenation (RF-lifting technology). The machine has robust wheels, which allows it's easy to move.

Operation principle

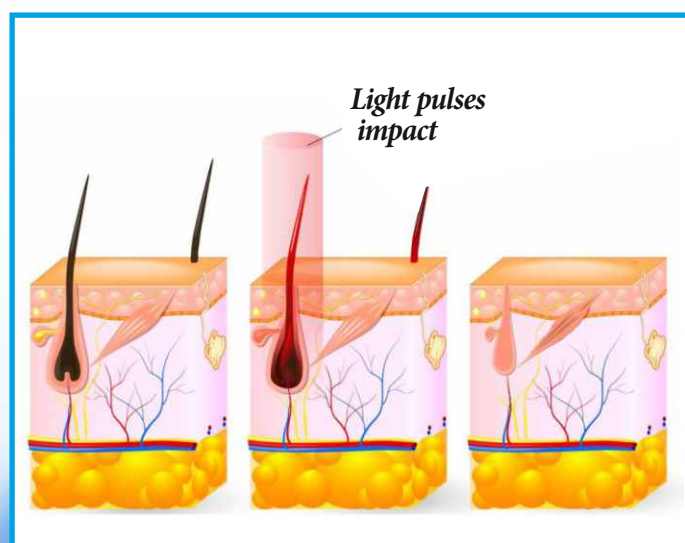

before melanin of hair shaft absorbs light and heats up complete destruction of the hair follicle

Advantages:

- Complete hair removal for long term
- Large spot size: 15 * 50 * 10 mm and 40 mm
- Pain-free
- Removes all hair types
- Suitable for all skin types
- Intelligent control regulating by innovative software

Applications:

- Hair Removal: Unwanted Hair (armpits, feet, lips, etc.)
- Skin rejuvenation using RF- technology lifting: pore reduction, tummy skin, improve skin elasticity

Parameters:

- Wavelength: HR: 690-950 nm, SR: 560-950 nm
- Energy density: HR: 1-10 J / cm², SR: 1-10 J / cm²
- Frequency: HR 1-8Gts, SR 1-2Gts
- Power: 3000 W
- Power Supply: pure sapphire
- spot size 10 * 40mm for the SR, 15 * 50 mm for HR
- Timer: 1s / 3s / 30s for HR, 1c / 3c / 10c for the SR
- Cooling system: integrated cooling system of the 4th degree
- Cooling system: water cooling + semiconductor cooling + air cooling
- LCD Screen: 10.4" touch screen



Body correction beauty machine PX-7000



Body correction beauty machine PX-7000+



Radiofrequency lifting with vacuum and cavitation PX-7000 beauty machine is an innovative system of rejuvenation on RF technology basis. Radiofrequency lifting beauty machine PX-7000 provides effective facial lifting and body shaping including RF correction and lifting, non-invasive ultrasonic liposuction and elimination of cellulite. The innovative technical achievements together with unique design allows the most effective improving of different body contours problems, modelling of body sculpture on the permanent basis comparable to plastic surgery. It is possible to use the machine functions as a complex and separately.



Parameters:

- nozzles: 1 cavitation nozzle, 1 vacuum nozzle with RF, 3 RF nozzles of different shape for different areas.
- LCD 7" screen
- Power: 100 W
- Cavitation frequency: 40 kHz
- RF frequency: 3-5 MHz
- Modes: pulse and continuous



Functions:

- cavitation 40K
- vacuum bipolar nozzle
- tripolar RF- tripolar nozzle
- body RF-multipolar nozzle
- face RF-multipolar nozzle

* PX-7000+ — the advantage of the machine is

Pressure therapy beauty machine PR - 801



Pressure therapy beauty

Air pressure body slimming beauty machines PR-801, PR-901 combine widely used body correction treatments such as:

- pressure therapy
- miostimulation
- infrared heating

Indications:

- activation of metabolic process
- acceleration of venous blood flow
- rehabilitation period after trauma
- edema
- varicosity
- elimination of toxins
- cellulitis
- overweight
- intercellular fluid circulation improvement

Functions:

- pressure therapy
- infrared heating
- muscle stimulation

Set of delivery:

- 20 air ports
- Body air bags -7 pcs.
- Blanket
- 5 bandages for muscle stimulation
- 10 electrode pads for muscle stimulation
- Heating of every part of suit up to 75° C




Pressure therapy beauty machine PR-2000



Parameters:

- pressure therapy machine with 48 air port
- intelligent interface
- 5 programmed modes
- mode of separate switching on / off sections
- individual settings of pressure intensity and time limit
- adjustable air supply: 1-6 sec.
- compressor endurance up to 18 hours per day
- working part of suit is displayed
- knee-high boots fixed on waiste
- long-sleeved suit fixed on breast

Pressure therapy beauty machine PR-701

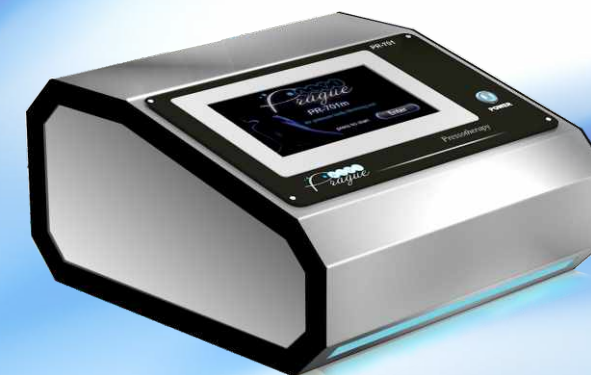

Air pressure body slimming beauty machine PR-701m is a professional device that has 16- channel suit, separated leg sectors with 4 ports on every leg, with 2 ports on every arm and 2 ports on the back area (small back and belly)

This machine is very lightweight and portable. It facilitates the transportation and storage of the beauty machine.



* Touch screen LCD 7 "

Assemble your own multifunctiona beauty machine!



*Beauty technology has
the name - it's **ALVI** Prague*

